CATTANI COMPRESSORS


AC400 5-6 CHAIRS

- Clean, Dry, Oil-Free Air
- 7-year Parts and Labour Warranty
- Outstanding Performance and Reliability
- Sound Reduction Options

Reliable: Robust Cattani compressors keep practices running, delivering clean, dry, oil-free air with such a level of reliability that we give it a 7-year warranty. The AC400, for 5-6 chairs, is easy to install and low maintenance.

Quality Protection: The tanks in all Cattani compressors have a patented silver ion coating that delivers highly effective protection against microbes such as bacteria and mould. Permanent antimicrobial protection is assured. Internal components are also specially coated to reduce mechanical wear and tear.

Sound Reduction: If the compressor is in close proximity to treatment and staff rooms, Cattani offers the AC400Q, which significantly reduces noise output.

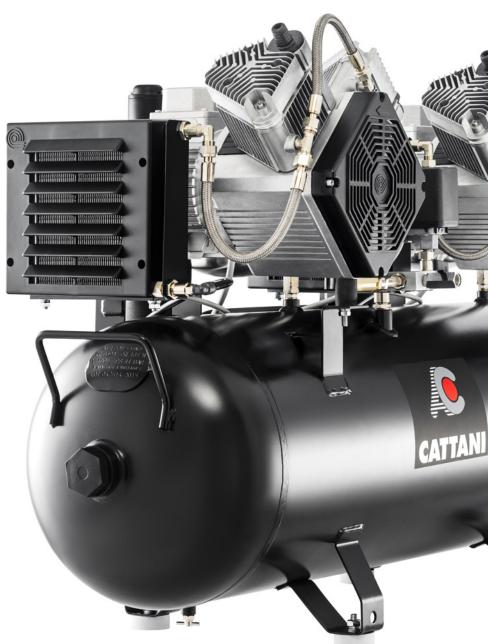



PERFORMANCE AND RELIABILITY

Cattani Australasia is a family company that stands for outstanding quality, performance and service. We have more than 30 years' experience in the field of dental air technology, and are proud to work with what we believe are the best compressors and suction systems in the business.

The confidence we have in our compressors' superior design, mechanical precision and robust components means we have no hesitation in backing them with an

industry-leading 7-year parts and labour warranty.



AC400 (THREE PHASE) & AC4001P (SINGLE PHASE)

NO. OF CHAIRS	5-6
POWER	3~ 400V 50Hz 7.4A 3kW 1~ 240V 50Hz 15.4A 2.4kW
HEIGHT MM	745
WIDTH MM	1155
DEPTH MM	550
WEIGHT KG	112
SOUND OUTPUT dB(A)	73
MAX AIR FLOW AT 5 BAR NL/min	320
TANK SIZE (LITRES)	100
PRODUCT CODES	
AC400	(AU/NZ) K-AC4000, (IT) 13432
AC4001P	(AU/NZ) K-AC4001P, (IT) 13422

Cattani compressors, which come in sizes to suit every clinic, are manufactured to the strictest standards in Parma, Italy. Each unit is rigorously tested before delivery, allowing the clinician to have the utmost confidence in its performance and reliability.

SOUND REDUCTION OPTIONS

In some clinics, plant room equipment noise can be an issue, for example where the plant room is close to treatment and staff rooms. In these situations, we have sound reduction options that can dramatically cut noise.



AC400Q (THREE PHASE) & AC400Q1P (SINGLE PHASE)

NO. OF CHAIRS	5-6
POWER	3~ 400V 50Hz 7.4A 3kW 1~ 230V 50Hz 15.4A 2.4kW
HEIGHT MM	890
WIDTH MM	1160
DEPTH MM	650
WEIGHT KG	135
SOUND OUTPUT dB(A)	63.8
MAX AIR FLOW AT 5 BAR NL/min	320
TANK SIZE (LITRES)	100
PRODUCT CODES AC400Q	(AU/NZ) K-AC400Q, (IT) 13450



VENTILATION KITS

Recommendations - 100 m3/h for each kW of equipment installed. 250 & 500 kits include fan. ducting and vent cover. 800 includes fan only.

Product codes:

250 - L-TD2500-KIT 500 - L-TD5000-150 800 - L-TD8000



HEPA H14 FILTER

For medical quality air. Non-autoclavable Two required.

Product code: K-168429



(AU/NZ) K-AC400Q1P, (IT) 13444

DRYER PRE FILTER FOR COMPRESSOR (FOR HUMID CLIMATES)

A dryer pre filter is recommended for humid climates. This additional moisture trap is fitted prior to the drying column. The pre filter has proven extremely effective and reliable in increasing the efficiency of the drying system by removing bulk condensation from the air prior to entering the drying column. Can be retro-fitted. Two required.

Product code: K-E10000

AC400Q1P

Australia

+61 3 9484 1120 cattani@cattani.com.au www.cattani.com.au

280 Dundas Street Thornbury VIC 3071

New Zealand

0800 68 22 88 cattani@cattani.co.nz www.cattani.co.nz

