

CATTANI AFTER CARE EQUIPMENT SERVICING SCHEDULE

CONTENTS

04

AFTER CARE PROGRAMME



- 5 How to promote to your customers
- 6 8 Ways to promote the After Care Programme to your customers

09

COMPRESSORS ANNUAL SERVICE



- 11 AC 100 - AC 200
- 19 AC 300 - AC 400
- 31 AC 600 - AC 900
- 37 AC 1200 - 1800
- 41 AC 2400
- 43 Block Jet 9 x 16 Cyl.

44

CAD CAM COMPRESSORS ANNUAL SERVICE



- 44 AC 210 - AC 310
- 48 AC 410 - AC 610
- 54 AC 910

55

SMART SERIES SUCTION SYSTEMS ANNUAL SERVICE



- 57 Micro SMART
- 59 Micro SMART Cube
- 61 Turbo SMART
- 65 Turbo SMART Cube

71

TANDEM SMART SERIES SUCTION SYSTEMS ANNUAL SERVICE



- 69 Tandem Micro SMART Range
- 71 Tandem Turbo SMART Range
- 75 Multi SMART T3 Range
- 77 Multi SMART T4 Range
- 79 Multi SMART T5 Range

81

MAXI SMART SERIES
SUCTION SYSTEMS
ANNUAL SERVICE



81 Maxi Smart

83

BLOCK JET LARGE
VACUUM SYSTEM
ANNUAL SERVICE



83 Block Jet Large Central
Vacuum System

85

SUCTION SYSTEMS
ANNUAL SERVICE



85 Aspi Jet Range
89 Turbo Jet Range

91

DRY SYSTEMS
ANNUAL SERVICE



81 Dry Suction Range



AFTER CARE PROGRAMME

CATTANI AFTER CARE – IT'S A SMART MOVE

How to promote the After Care Programme to your customers

Cattani After Care is a service programme with genuine benefits for dentists. This document lists the key advantages and in simple language and suggests different ways to express those advantages.

You won't need to use all of the arguments presented here, so choose the ones that suit you, the dentist and the practice situation.

And don't hesitate to speak to your friends at Cattani if you need any more information or advice regarding the programme.

8 Ways to promote the After Care Programme to your customers

1. Your annual check-up is due

Preventive maintenance is a proven concept – even dentists do it

Dentists encourage their patients to have regular check-ups and oral hygiene visits because it helps them keep their natural teeth and delays major work for as long as possible. Preventive machinery maintenance like the Cattani After Care Programme does the same thing. It's like a regular health check for the dentist's plant room equipment.

2. Scheduled servicing keeps a motor running

All engines need regular maintenance – it's a mechanical fact of life

A Cattani machine is more like a car than any other piece of equipment in a dental practice, and like a car it's unrealistic to think it will run indefinitely without some TLC. No one buys a Ferrari and runs it into the ground – same with Cattani. It will perform longer, and much more reliably, if it is looked after on a regular basis with the Cattani After Care Programme.

3. It's basic risk management

Minimise the risk of breakdowns and their unpleasant consequences

Regular scheduled servicing of Cattani suction units and compressors minimises the risk of them breaking down. This idea can also be framed as “prevent now or pay later” or even as a kind of insurance. Think about reducing the chance of:

- Premature equipment failure
- Unscheduled down time – when is it ever be convenient?
- Expensive, unscheduled technician call-outs
- Significant machinery repair and replacement costs
- Embarrassment in front of patients – and the cost and inconvenience of rescheduling, maybe even losing a patient

4. Cattani After Care means no surprises

The Cattani After Care Programme is transparent in concept and cost

Cattani After Care is more than just regular scheduled servicing. Think of receiving an invoice after the annual car service, not knowing whether it's going to cost, say, \$300 or \$1200. That's another reason why Cattani After Care makes such good sense – you know what each service costs, so no bill shock. And every service in the programme has its own unique, pre-packaged kit with all the essential parts.

5. Cattani time is productive time

Cattani After Care keeps practice disruption to a minimum

This is a completely proactive service programme, managed by the technician. In many cases the technician can service Cattani equipment at the same time as other surgery equipment like autoclaves, chairs, etc, with visits coordinated to minimise costs and inconvenience for both parties.

6. Cattani After Care is in expert hands

No other specialist plant room equipment supplier matches Cattani's level of service

Cattani's intention is to provide the best level of service to end users (and service training courses for technicians), and dentists can be assured that only technicians fully trained in servicing Cattani products are invited to be part of the After Care Programme.

7. Everyone wins

Quality, documented service every time with Cattani After Care

Cattani After Care is a fully laid out programme with dates and terms and conditions, all in line with factory-recommended service intervals. After Care comes with a simple log book that is kept up to date by the technician. This log book leaves no doubt about where things are at with a particular machine.

8. Quality equipment, quality service, always

Optimal performance and professionalism is good for everyone

No one can afford to be casual about their professional reputations. Cattani After Care is an easy way to protect and even enhance the reputations of the dentist, the technician and Cattani itself. We all take a hit if something goes wrong, so let's do everything we can to prevent that, and in fact go the extra mile. Grow our reputations and we grow our businesses.

AC SERIES COMPRESSORS



AC SERIES COMPRESSOR ANNUAL SERVICING SCHEDULE

Cattani Compressors are serviced annually

- 1x Service Kit required per head unit
- Kits are available for head units with 1,2, 3 & 6 cylinder(s)
- All Service Kits available with a standard filter or HEPA filter for the drying column

AC COMPRESSOR (~1)

Annual Servicing Schedule

Year 1 & Year 2 - Minor Service Kit

- Air inlet filters (1 filter per cylinder), drying column filter and o-ring seal
- HEPA filter option for drying column is available.

Year 3 - Intermediate Service Kit

- Air inlet filters (1 filter per cylinder), drying column filter, o-ring seal and a motor start capacitor
- HEPA filter option for drying column is available.

AC COMPRESSOR (~3)

Annual Servicing Schedule

Every Year - Minor Service Kit

- Air inlet filters (1x filter per cylinder), drying column filter and o-ring seal
- HEPA filter option for drying column is available.

ANNUAL SERVICING SCHEDULE

AC100 Compressor (~1)

(with Standard Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK1013 Inlet filters & O-ring	
2	(x1) K-SK1013 Inlet filters & O-ring	
3		(x1) K-SK1010 Capacitor, Inlet filters & O-ring
4	(x1) K-SK1013 Inlet filters & O-ring	
5	(x1) K-SK1013 Inlet filters & O-ring	
6		(x1) K-SK1010 Capacitor, Inlet filters & O-ring
7	(x1) K-SK1013 Inlet filters & O-ring	
8	(x1) K-SK1013 Inlet filters & O-ring	
9		(x1) K-SK1010 Capacitor, Inlet filters & O-ring
10	(x1) K-SK1013 Inlet filters & O-ring	

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC100 Compressor (~1)

(with HEPA Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK1013B Inlet filters & O-ring	
2	(x1) K-SK1013B Inlet filters & O-ring	
3		(x1) K-SK101B Capacitor, Inlet filters & O-ring
4	(x1) K-SK1013B Inlet filters & O-ring	
5	(x1) K-SK1013B Inlet filters & O-ring	
6		(x1) K-SK101B Capacitor, Inlet filters & O-ring
7	(x1) K-SK1013B Inlet filters & O-ring	
8	(x1) K-SK1013B Inlet filters & O-ring	
9		(x1) K-SK101B Capacitor, Inlet filters & O-ring
10	(x1) K-SK1013B Inlet filters & O-ring	

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC100Q Compressor (~1)

(with Standard Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK1013 Inlet filters & O-ring	
2	(x1) K-SK1013 Inlet filters & O-ring	
3		(x1) K-SK1010 Capacitor, Inlet filters & O-ring
4	(x1) K-SK1013 Inlet filters & O-ring	
5	(x1) K-SK1013 Inlet filters & O-ring	
6		(x1) K-SK1010 Capacitor, Inlet filters & O-ring
7	(x1) K-SK1013 Inlet filters & O-ring	
8	(x1) K-SK1013 Inlet filters & O-ring	
9		(x1) K-SK1010 Capacitor, Inlet filters & O-ring
10	(x1) K-SK1013 Inlet filters & O-ring	

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC100Q Compressor (~1)

(with HEPA Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK1013B Inlet filters & O-ring	
2	(x1) K-SK1013B Inlet filters & O-ring	
3		(x1) K-SK101B Capacitor, Inlet filters & O-ring
4	(x1) K-SK1013B Inlet filters & O-ring	
5	(x1) K-SK1013B Inlet filters & O-ring	
6		(x1) K-SK101B Capacitor, Inlet filters & O-ring
7	(x1) K-SK1013B Inlet filters & O-ring	
8	(x1) K-SK1013B Inlet filters & O-ring	
9		(x1) K-SK101B Capacitor, Inlet filters & O-ring
10	(x1) K-SK1013B Inlet filters & O-ring	

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC200 Compressor (~1)

(with Standard Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK2013 Inlet filters & O-ring	
2	(x1) K-SK2013 Inlet filters & O-ring	
3		(x1) K-SK2010 Capacitor, Inlet filters & O-ring
4	(x1) K-SK2013 Inlet filters & O-ring	
5	(x1) K-SK2013 Inlet filters & O-ring	
6		(x1) K-SK2010 Capacitor, Inlet filters & O-ring
7	(x1) K-SK2013 Inlet filters & O-ring	
8	(x1) K-SK2013 Inlet filters & O-ring	
9		(x1) K-SK2010 Capacitor, Inlet filters & O-ring
10	(x1) K-SK2013 Inlet filters & O-ring	

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC200 Compressor (~1)

(with HEPA Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK2013B Inlet filters & O-ring	
2	(x1) K-SK2013B Inlet filters & O-ring	
3		(x1) K-SK201B Capacitor, Inlet filters & O-ring
4	(x1) K-SK2013B Inlet filters & O-ring	
5	(x1) K-SK2013B Inlet filters & O-ring	
6		(x1) K-SK201B Capacitor, Inlet filters & O-ring
7	(x1) K-SK2013B Inlet filters & O-ring	
8	(x1) K-SK2013B Inlet filters & O-ring	
9		(x1) K-SK201B Capacitor, Inlet filters & O-ring
10	(x1) K-SK2013B Inlet filters & O-ring	

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC200Q Compressor (~1)

(with Standard Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK2013 Inlet filters & O-ring	
2	(x1) K-SK2013 Inlet filters & O-ring	
3		(x1) K-SK2010 Capacitor, Inlet filters & O-ring
4	(x1) K-SK2013 Inlet filters & O-ring	
5	(x1) K-SK2013 Inlet filters & O-ring	
6		(x1) K-SK2010 Capacitor, Inlet filters & O-ring
7	(x1) K-SK2013 Inlet filters & O-ring	
8	(x1) K-SK2013 Inlet filters & O-ring	
9		(x1) K-SK2010 Capacitor, Inlet filters & O-ring
10	(x1) K-SK2013 Inlet filters & O-ring	

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC200Q Compressor (~1)

(with HEPA Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK2013B Inlet filters & O-ring	
2	(x1) K-SK2013B Inlet filters & O-ring	
3		(x1) K-SK201B Capacitor, Inlet filters & O-ring
4	(x1) K-SK2013B Inlet filters & O-ring	
5	(x1) K-SK2013B Inlet filters & O-ring	
6		(x1) K-SK201B Capacitor, Inlet filters & O-ring
7	(x1) K-SK2013B Inlet filters & O-ring	
8	(x1) K-SK2013B Inlet filters & O-ring	
9		(x1) K-SK201B Capacitor, Inlet filters & O-ring
10	(x1) K-SK2013B Inlet filters & O-ring	

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC300 Compressor (~1)

(with Standard Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK3013 Inlet filters & O-ring	
2	(x1) K-SK3013 Inlet filters & O-ring	
3		(x1) K-SK3010 Capacitor, Inlet filters & O-ring
4	(x1) K-SK3013 Inlet filters & O-ring	
5	(x1) K-SK3013 Inlet filters & O-ring	
6		(x1) K-SK3010 Capacitor, Inlet filters & O-ring
7	(x1) K-SK3013 Inlet filters & O-ring	
8	(x1) K-SK3013 Inlet filters & O-ring	
9		(x1) K-SK3010 Capacitor, Inlet filters & O-ring
10	(x1) K-SK3013 Inlet filters & O-ring	

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC300 Compressor (~1)

(with HEPA Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK3013B Inlet filters & O-ring	
2	(x1) K-SK3013B Inlet filters & O-ring	
3		(x1) K-SK301B Capacitor, Inlet filters & O-ring
4	(x1) K-SK3013B Inlet filters & O-ring	
5	(x1) K-SK3013B Inlet filters & O-ring	
6		(x1) K-SK301B Capacitor, Inlet filters & O-ring
7	(x1) K-SK3013B Inlet filters & O-ring	
8	(x1) K-SK3013B Inlet filters & O-ring	
9		(x1) K-SK301B Capacitor, Inlet filters & O-ring
10	(x1) K-SK3013B Inlet filters & O-ring	

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC300Q Compressor (~1)

(with Standard Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK3013 Inlet filters & O-ring	
2	(x1) K-SK3013 Inlet filters & O-ring	
3		(x1) K-SK3010 Capacitor, Inlet filters & O-ring
4	(x1) K-SK3013 Inlet filters & O-ring	
5	(x1) K-SK3013 Inlet filters & O-ring	
6		(x1) K-SK3010 Capacitor, Inlet filters & O-ring
7	(x1) K-SK3013 Inlet filters & O-ring	
8	(x1) K-SK3013 Inlet filters & O-ring	
9		(x1) K-SK3010 Capacitor, Inlet filters & O-ring
10	(x1) K-SK3013 Inlet filters & O-ring	

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC300Q Compressor (~1)

(with HEPA Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK3013B Inlet filters & O-ring	
2	(x1) K-SK3013B Inlet filters & O-ring	
3		(x1) K-SK301B Capacitor, Inlet filters & O-ring
4	(x1) K-SK3013B Inlet filters & O-ring	
5	(x1) K-SK3013B Inlet filters & O-ring	
6		(x1) K-SK301B Capacitor, Inlet filters & O-ring
7	(x1) K-SK3013B Inlet filters & O-ring	
8	(x1) K-SK3013B Inlet filters & O-ring	
9		(x1) K-SK301B Capacitor, Inlet filters & O-ring
10	(x1) K-SK3013B Inlet filters & O-ring	

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC400P Compressor (~1)

(with Standard Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x2) K-SK2013 Inlet filters & O-ring	
2	(x2) K-SK2013 Inlet filters & O-ring	
3		(x2) K-SK2010 Capacitor, Inlet filters & O-ring
4	(x2) K-SK2013 Inlet filters & O-ring	
5	(x2) K-SK2013 Inlet filters & O-ring	
6		(x2) K-SK2010 Capacitor, Inlet filters & O-ring
7	(x2) K-SK2013 Inlet filters & O-ring	
8	(x2) K-SK2013 Inlet filters & O-ring	
9		(x2) K-SK2010 Capacitor, Inlet filters & O-ring
10	(x2) K-SK2013 Inlet filters & O-ring	

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC4001P Compressor (~1)

(with HEPA Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x2) K-SK2013B Inlet filters & O-ring	
2	(x2) K-SK2013B Inlet filters & O-ring	
3		(x2) K-SK201B Capacitor, Inlet filters & O-ring
4	(x2) K-SK2013B Inlet filters & O-ring	
5	(x2) K-SK2013B Inlet filters & O-ring	
6		(x2) K-SK201B Capacitor, Inlet filters & O-ring
7	(x2) K-SK2013B Inlet filters & O-ring	
8	(x2) K-SK2013B Inlet filters & O-ring	
9		(x2) K-SK201B Capacitor, Inlet filters & O-ring
10	(x2) K-SK2013B Inlet filters & O-ring	

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC400Q1P Compressor (~1)

(with Standard Filter for drying column)



* Servicing schedule also applies to “E” version of this compressor

Year	Minor Service Kits	Intermediate Service Kits
1	(x2) K-SK2013 Inlet filters & O-ring	
2	(x2) K-SK2013 Inlet filters & O-ring	
3		(x2) K-SK2010 Capacitor, Inlet filters & O-ring
4	(x2) K-SK2013 Inlet filters & O-ring	
5	(x2) K-SK2013 Inlet filters & O-ring	
6		(x2) K-SK2010 Capacitor, Inlet filters & O-ring
7	(x2) K-SK2013 Inlet filters & O-ring	
8	(x2) K-SK2013 Inlet filters & O-ring	
9		(x2) K-SK2010 Capacitor, Inlet filters & O-ring
10	(x2) K-SK2013 Inlet filters & O-ring	

ANNUAL SERVICING SCHEDULE

AC400Q1P Compressor (~1)

(with HEPA Filter for drying column)



* Servicing schedule also applies to “E” version of this compressor

Year	Minor Service Kits	Intermediate Service Kits
1	(x2) K-SK2013B Inlet filters & O-ring	
2	(x2) K-SK2013B Inlet filters & O-ring	
3		(x2) K-SK201B Capacitor, Inlet filters & O-ring
4	(x2) K-SK2013B Inlet filters & O-ring	
5	(x2) K-SK2013B Inlet filters & O-ring	
6		(x2) K-SK201B Capacitor, Inlet filters & O-ring
7	(x2) K-SK2013B Inlet filters & O-ring	
8	(x2) K-SK2013B Inlet filters & O-ring	
9		(x2) K-SK201B Capacitor, Inlet filters & O-ring
10	(x2) K-SK2013B Inlet filters & O-ring	

ANNUAL SERVICING SCHEDULE

AC400 Compressor (~3)

(with Standard Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK2013 Inlet filters & O-ring
2	(x2) K-SK2013 Inlet filters & O-ring
3	(x2) K-SK2013 Inlet filters & O-ring
4	(x2) K-SK2013 Inlet filters & O-ring
5	(x2) K-SK2013 Inlet filters & O-ring
6	(x2) K-SK2013 Inlet filters & O-ring
7	(x2) K-SK2013 Inlet filters & O-ring
8	(x2) K-SK2013 Inlet filters & O-ring
9	(x2) K-SK2013 Inlet filters & O-ring
10	(x2) K-SK2013 Inlet filters & O-ring

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC400 Compressor (~3)

(with HEPA Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK2013B Inlet filters & O-ring
2	(x2) K-SK2013B Inlet filters & O-ring
3	(x2) K-SK2013B Inlet filters & O-ring
4	(x2) K-SK2013B Inlet filters & O-ring
5	(x2) K-SK2013B Inlet filters & O-ring
6	(x2) K-SK2013B Inlet filters & O-ring
7	(x2) K-SK2013B Inlet filters & O-ring
8	(x2) K-SK2013B Inlet filters & O-ring
9	(x2) K-SK2013B Inlet filters & O-ring
10	(x2) K-SK2013B Inlet filters & O-ring

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC400Q Compressor (~3)

(with Standard Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK2013 Inlet filters & O-ring
2	(x2) K-SK2013 Inlet filters & O-ring
3	(x2) K-SK2013 Inlet filters & O-ring
4	(x2) K-SK2013 Inlet filters & O-ring
5	(x2) K-SK2013 Inlet filters & O-ring
6	(x2) K-SK2013 Inlet filters & O-ring
7	(x2) K-SK2013 Inlet filters & O-ring
8	(x2) K-SK2013 Inlet filters & O-ring
9	(x2) K-SK2013 Inlet filters & O-ring
10	(x2) K-SK2013 Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC400Q Compressor (~3)

(with HEPA Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK2013B Inlet filters & O-ring
2	(x2) K-SK2013B Inlet filters & O-ring
3	(x2) K-SK2013B Inlet filters & O-ring
4	(x2) K-SK2013B Inlet filters & O-ring
5	(x2) K-SK2013B Inlet filters & O-ring
6	(x2) K-SK2013B Inlet filters & O-ring
7	(x2) K-SK2013B Inlet filters & O-ring
8	(x2) K-SK2013B Inlet filters & O-ring
9	(x2) K-SK2013B Inlet filters & O-ring
10	(x2) K-SK2013B Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC600 Compressor (~3)

(with Standard Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK3013 Inlet filters & O-ring
2	(x2) K-SK3013 Inlet filters & O-ring
3	(x2) K-SK3013 Inlet filters & O-ring
4	(x2) K-SK3013 Inlet filters & O-ring
5	(x2) K-SK3013 Inlet filters & O-ring
6	(x2) K-SK3013 Inlet filters & O-ring
7	(x2) K-SK3013 Inlet filters & O-ring
8	(x2) K-SK3013 Inlet filters & O-ring
9	(x2) K-SK3013 Inlet filters & O-ring
10	(x2) K-SK3013 Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC600 Compressor (~3)

(with HEPA Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK3013B Inlet filters & O-ring
2	(x2) K-SK3013B Inlet filters & O-ring
3	(x2) K-SK3013B Inlet filters & O-ring
4	(x2) K-SK3013B Inlet filters & O-ring
5	(x2) K-SK3013B Inlet filters & O-ring
6	(x2) K-SK3013B Inlet filters & O-ring
7	(x2) K-SK3013B Inlet filters & O-ring
8	(x2) K-SK3013B Inlet filters & O-ring
9	(x2) K-SK3013B Inlet filters & O-ring
10	(x2) K-SK3013B Inlet filters & O-ring

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC600Q Compressor (~3)

(with Standard Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK3013 Inlet filters & O-ring
2	(x2) K-SK3013 Inlet filters & O-ring
3	(x2) K-SK3013 Inlet filters & O-ring
4	(x2) K-SK3013 Inlet filters & O-ring
5	(x2) K-SK3013 Inlet filters & O-ring
6	(x2) K-SK3013 Inlet filters & O-ring
7	(x2) K-SK3013 Inlet filters & O-ring
8	(x2) K-SK3013 Inlet filters & O-ring
9	(x2) K-SK3013 Inlet filters & O-ring
10	(x2) K-SK3013 Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC600Q Compressor (~3)

(with HEPA Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK3013B Inlet filters & O-ring
2	(x2) K-SK3013B Inlet filters & O-ring
3	(x2) K-SK3013B Inlet filters & O-ring
4	(x2) K-SK3013B Inlet filters & O-ring
5	(x2) K-SK3013B Inlet filters & O-ring
6	(x2) K-SK3013B Inlet filters & O-ring
7	(x2) K-SK3013B Inlet filters & O-ring
8	(x2) K-SK3013B Inlet filters & O-ring
9	(x2) K-SK3013B Inlet filters & O-ring
10	(x2) K-SK3013B Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC900 Compressor (~3)

(with Standard Filter for drying column)



Year	Minor Service Kits
1	(x3) K-SK3013 Inlet filters & O-ring
2	(x3) K-SK3013 Inlet filters & O-ring
3	(x3) K-SK3013 Inlet filters & O-ring
4	(x3) K-SK3013 Inlet filters & O-ring
5	(x3) K-SK3013 Inlet filters & O-ring
6	(x3) K-SK3013 Inlet filters & O-ring
7	(x3) K-SK3013 Inlet filters & O-ring
8	(x3) K-SK3013 Inlet filters & O-ring
9	(x3) K-SK3013 Inlet filters & O-ring
10	(x3) K-SK3013 Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC900 Compressor (~3)

(with HEPA Filter for drying column)



Year	Minor Service Kits
1	(x3) K-SK3013B Inlet filters & O-ring
2	(x3) K-SK3013B Inlet filters & O-ring
3	(x3) K-SK3013B Inlet filters & O-ring
4	(x3) K-SK3013B Inlet filters & O-ring
5	(x3) K-SK3013B Inlet filters & O-ring
6	(x3) K-SK3013B Inlet filters & O-ring
7	(x3) K-SK3013B Inlet filters & O-ring
8	(x3) K-SK3013B Inlet filters & O-ring
9	(x3) K-SK3013B Inlet filters & O-ring
10	(x3) K-SK3013B Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC1200 Compressor (~3)

(with Standard Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK0106-11 Inlet filters & O-ring
2	(x2) K-SK0106-11 Inlet filters & O-ring
3	(x2) K-SK0106-11 Inlet filters & O-ring
4	(x2) K-SK0106-11 Inlet filters & O-ring
5	(x2) K-SK0106-11 Inlet filters & O-ring
6	(x2) K-SK0106-11 Inlet filters & O-ring
7	(x2) K-SK0106-11 Inlet filters & O-ring
8	(x2) K-SK0106-11 Inlet filters & O-ring
9	(x2) K-SK0106-11 Inlet filters & O-ring
10	(x2) K-SK0106-11 Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC1200 Compressor (~3)

(with HEPA Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK0106-11B Inlet filters & O-ring
2	(x2) K-SK0106-11B Inlet filters & O-ring
3	(x2) K-SK0106-11B Inlet filters & O-ring
4	(x2) K-SK0106-11B Inlet filters & O-ring
5	(x2) K-SK0106-11B Inlet filters & O-ring
6	(x2) K-SK0106-11B Inlet filters & O-ring
7	(x2) K-SK0106-11B Inlet filters & O-ring
8	(x2) K-SK0106-11B Inlet filters & O-ring
9	(x2) K-SK0106-11B Inlet filters & O-ring
10	(x2) K-SK0106-11B Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC1800 Compressor (~3)

(with Standard Filter for drying column)



Year	Minor Service Kits
1	(x3) K-SK0106-11 Inlet filters & O-ring
2	(x3) K-SK0106-11 Inlet filters & O-ring
3	(x3) K-SK0106-11 Inlet filters & O-ring
4	(x3) K-SK0106-11 Inlet filters & O-ring
5	(x3) K-SK0106-11 Inlet filters & O-ring
6	(x3) K-SK0106-11 Inlet filters & O-ring
7	(x3) K-SK0106-11 Inlet filters & O-ring
8	(x3) K-SK0106-11 Inlet filters & O-ring
9	(x3) K-SK0106-11 Inlet filters & O-ring
10	(x3) K-SK0106-11 Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC1800 Compressor (~3)

(with HEPA Filter for drying column)



Year	Minor Service Kits
1	(x3) K-SK0106-11B Inlet filters & O-ring
2	(x3) K-SK0106-11B Inlet filters & O-ring
3	(x3) K-SK0106-11B Inlet filters & O-ring
4	(x3) K-SK0106-11B Inlet filters & O-ring
5	(x3) K-SK0106-11B Inlet filters & O-ring
6	(x3) K-SK0106-11B Inlet filters & O-ring
7	(x3) K-SK0106-11B Inlet filters & O-ring
8	(x3) K-SK0106-11B Inlet filters & O-ring
9	(x3) K-SK0106-11B Inlet filters & O-ring
10	(x3) K-SK0106-11B Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC2400 Compressor (~3)

(with Standard Filtration)



Year	Minor Service Kits
1	(x24) K-168563 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
2	(x24) K-168563 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
3	(x24) K-168563 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
4	(x24) K-168563 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
5	(x24) K-168563 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)

Year	Minor Service Kits
6	(x24) K-168563 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
7	(x24) K-168563 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
8	(x24) K-168563 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
9	(x24) K-168563 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
10	(x24) K-168563 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)

ANNUAL SERVICING SCHEDULE

AC2400 Compressor (~3)

(with HEPA H14 Filtration)



Year	Minor Service Kits
1	(x24) K-168429 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
2	(x24) K-168429 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
3	(x24) K-168429 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
4	(x24) K-168429 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
5	(x24) K-168429 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)

Year	Minor Service Kits
6	(x24) K-168429 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
7	(x24) K-168429 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
8	(x24) K-168429 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
9	(x24) K-168429 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)
10	(x24) K-168429 (Inlet Filters) (x1) C-ZV2020 (3 Micron Prefilter) (x1) C-ZX2020 (0.01 Micron Filter)

ANNUAL SERVICING SCHEDULE

Blockjet (9x) 6 Cylinder

(with HEPA H14 Filtration)



Year	Minor Service Kits
1	(x54) K-168429 (Inlet Filters) (x1) F-211771 (3 Micron Pre-Filter) (x1) F-211773 (0.01 Micron Post-Filter)
2	(x54) K-168429 (Inlet Filters) (x1) F-211771 (3 Micron Pre-Filter) (x1) F-211773 (0.01 Micron Post-Filter)
3	(x54) K-168429 (Inlet Filters) (x1) F-211771 (3 Micron Pre-Filter) (x1) F-211773 (0.01 Micron Post-Filter)
4	(x54) K-168429 (Inlet Filters) (x1) F-211771 (3 Micron Pre-Filter) (x1) F-211773 (0.01 Micron Post-Filter)
5	(x54) K-168429 (Inlet Filters) (x1) F-211771 (3 Micron Pre-Filter) (x1) F-211773 (0.01 Micron Post-Filter)

Year	Minor Service Kits
6	(x54) K-168429 (Inlet Filters) (x1) F-211771 (3 Micron Pre-Filter) (x1) F-211773 (0.01 Micron Post-Filter)
7	(x54) K-168429 (Inlet Filters) (x1) F-211771 (3 Micron Pre-Filter) (x1) F-211773 (0.01 Micron Post-Filter)
8	(x54) K-168429 (Inlet Filters) (x1) F-211771 (3 Micron Pre-Filter) (x1) F-211773 (0.01 Micron Post-Filter)
9	(x54) K-168429 (Inlet Filters) (x1) F-211771 (3 Micron Pre-Filter) (x1) F-211773 (0.01 Micron Post-Filter)
10	(x54) K-168429 (Inlet Filters) (x1) F-211771 (3 Micron Pre-Filter) (x1) F-211773 (0.01 Micron Post-Filter)

ANNUAL SERVICING SCHEDULE

AC210 Compressor (~1)

(with HEPA Filter for drying column)



Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK2013B Inlet filters & O-ring	
2	(x1) K-SK2013B Inlet filters & O-ring	
3		(x1) K-SK201B Capacitor, Inlet filters & O-ring
4	(x1) K-SK2013B Inlet filters & O-ring	
5	(x1) K-SK2013B Inlet filters & O-ring	
6		(x1) K-SK201B Capacitor, Inlet filters & O-ring
7	(x1) K-SK2013B Inlet filters & O-ring	
8	(x1) K-SK2013B Inlet filters & O-ring	
9		(x1) K-SK201B Capacitor, Inlet filters & O-ring
10	(x1) K-SK2013B Inlet filters & O-ring	

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC210Q Compressor (~1)

(with HEPA Filter for drying column)



* Servicing schedule also applies to “E” version of this compressor

Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK2013B Inlet filters & O-ring	
2	(x1) K-SK2013B Inlet filters & O-ring	
3		(x1) K-SK201B Capacitor, Inlet filters & O-ring
4	(x1) K-SK2013B Inlet filters & O-ring	
5	(x1) K-SK2013B Inlet filters & O-ring	
6		(x1) K-SK201B Capacitor, Inlet filters & O-ring
7	(x1) K-SK2013B Inlet filters & O-ring	
8	(x1) K-SK2013B Inlet filters & O-ring	
9		(x1) K-SK201B Capacitor, Inlet filters & O-ring
10	(x1) K-SK2013B Inlet filters & O-ring	

ANNUAL SERVICING SCHEDULE

AC310 Compressor (~1)

(with HEPA Filter for drying column)



* Servicing schedule also applies to "E" version of this compressor

Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK3013B Inlet filters & O-ring	
2	(x1) K-SK3013B Inlet filters & O-ring	
3		(x1) K-SK301B Capacitor, Inlet filters & O-ring
4	(x1) K-SK3013B Inlet filters & O-ring	
5	(x1) K-SK3013B Inlet filters & O-ring	
6		(x1) K-SK301B Capacitor, Inlet filters & O-ring
7	(x1) K-SK3013B Inlet filters & O-ring	
8	(x1) K-SK3013B Inlet filters & O-ring	
9		(x1) K-SK301B Capacitor, Inlet filters & O-ring
10	(x1) K-SK3013B Inlet filters & O-ring	

ANNUAL SERVICING SCHEDULE

AC310Q Compressor (~1)

(with HEPA Filter for drying column)



* Servicing schedule also applies to "E" version of this compressor

Year	Minor Service Kits	Intermediate Service Kits
1	(x1) K-SK3013B Inlet filters & O-ring	
2	(x1) K-SK3013B Inlet filters & O-ring	
3		(x1) K-SK301B Capacitor, Inlet filters & O-ring
4	(x1) K-SK3013B Inlet filters & O-ring	
5	(x1) K-SK3013B Inlet filters & O-ring	
6		(x1) K-SK301B Capacitor, Inlet filters & O-ring
7	(x1) K-SK3013B Inlet filters & O-ring	
8	(x1) K-SK3013B Inlet filters & O-ring	
9		(x1) K-SK301B Capacitor, Inlet filters & O-ring
10	(x1) K-SK3013B Inlet filters & O-ring	

ANNUAL SERVICING SCHEDULE

AC4101P Compressor (~1)



Year	Minor Service Kits	Intermediate Service Kits
1	(x2) K-SK2013B Inlet filters & O-ring	
2	(x2) K-SK2013B Inlet filters & O-ring	
3		(x2) K-SK201B Capacitor, Inlet filters & O-ring
4	(x2) K-SK2013B Inlet filters & O-ring	
5	(x2) K-SK2013B Inlet filters & O-ring	
6		(x2) K-SK201B Capacitor, Inlet filters & O-ring
7	(x2) K-SK2013B Inlet filters & O-ring	
8	(x2) K-SK2013B Inlet filters & O-ring	
9		(x2) K-SK201B Capacitor, Inlet filters & O-ring
10	(x2) K-SK2013B Inlet filters & O-ring	

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC410Q1P Compressor (~1)



* Servicing schedule also applies to “E” version of this compressor

Year	Minor Service Kits	Intermediate Service Kits
1	(x2) K-SK2013B Inlet filters & O-ring	
2	(x2) K-SK2013B Inlet filters & O-ring	
3		(x2) K-SK201B Capacitor, Inlet filters & O-ring
4	(x2) K-SK2013B Inlet filters & O-ring	
5	(x2) K-SK2013B Inlet filters & O-ring	
6		(x2) K-SK201B Capacitor, Inlet filters & O-ring
7	(x2) K-SK2013B Inlet filters & O-ring	
8	(x2) K-SK2013B Inlet filters & O-ring	
9		(x2) K-SK201B Capacitor, Inlet filters & O-ring
10	(x2) K-SK2013B Inlet filters & O-ring	

ANNUAL SERVICING SCHEDULE

AC410 Compressor (~3)

(with HEPA Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK2013B Inlet filters & O-ring
2	(x2) K-SK2013B Inlet filters & O-ring
3	(x2) K-SK2013B Inlet filters & O-ring
4	(x2) K-SK2013B Inlet filters & O-ring
5	(x2) K-SK2013B Inlet filters & O-ring
6	(x2) K-SK2013B Inlet filters & O-ring
7	(x2) K-SK2013B Inlet filters & O-ring
8	(x2) K-SK2013B Inlet filters & O-ring
9	(x2) K-SK2013B Inlet filters & O-ring
10	(x2) K-SK2013B Inlet filters & O-ring

* Servicing schedule also applies to "E" version of this compressor

ANNUAL SERVICING SCHEDULE

AC410Q Compressor (~3)

(with HEPA Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK2013B Inlet filters & O-ring
2	(x2) K-SK2013B Inlet filters & O-ring
3	(x2) K-SK2013B Inlet filters & O-ring
4	(x2) K-SK2013B Inlet filters & O-ring
5	(x2) K-SK2013B Inlet filters & O-ring
6	(x2) K-SK2013B Inlet filters & O-ring
7	(x2) K-SK2013B Inlet filters & O-ring
8	(x2) K-SK2013B Inlet filters & O-ring
9	(x2) K-SK2013B Inlet filters & O-ring
10	(x2) K-SK2013B Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC610 Compressor (~3)

(with HEPA Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK3013B Inlet filters & O-ring
2	(x2) K-SK3013B Inlet filters & O-ring
3	(x2) K-SK3013B Inlet filters & O-ring
4	(x2) K-SK3013B Inlet filters & O-ring
5	(x2) K-SK3013B Inlet filters & O-ring
6	(x2) K-SK3013B Inlet filters & O-ring
7	(x2) K-SK3013B Inlet filters & O-ring
8	(x2) K-SK3013B Inlet filters & O-ring
9	(x2) K-SK3013B Inlet filters & O-ring
10	(x2) K-SK3013B Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC610Q Compressor (~3)

(with HEPA Filter for drying column)



Year	Minor Service Kits
1	(x2) K-SK3013B Inlet filters & O-ring
2	(x2) K-SK3013B Inlet filters & O-ring
3	(x2) K-SK3013B Inlet filters & O-ring
4	(x2) K-SK3013B Inlet filters & O-ring
5	(x2) K-SK3013B Inlet filters & O-ring
6	(x2) K-SK3013B Inlet filters & O-ring
7	(x2) K-SK3013B Inlet filters & O-ring
8	(x2) K-SK3013B Inlet filters & O-ring
9	(x2) K-SK3013B Inlet filters & O-ring
10	(x2) K-SK3013B Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor

ANNUAL SERVICING SCHEDULE

AC910 Compressor (~3)

(with HEPA Filter for drying column)



Year	Minor Service Kits
1	(x3) K-SK3013B Inlet filters & O-ring
2	(x3) K-SK3013B Inlet filters & O-ring
3	(x3) K-SK3013B Inlet filters & O-ring
4	(x3) K-SK3013B Inlet filters & O-ring
5	(x3) K-SK3013B Inlet filters & O-ring
6	(x3) K-SK3013B Inlet filters & O-ring
7	(x3) K-SK3013B Inlet filters & O-ring
8	(x3) K-SK3013B Inlet filters & O-ring
9	(x3) K-SK3013B Inlet filters & O-ring
10	(x3) K-SK3013B Inlet filters & O-ring

* Servicing schedule also applies to “E” version of this compressor



SMART SERIES SUCTION



SMART SERIES SUCTION ANNUAL SERVICING SCHEDULE

Cattani SMART Suction Units are serviced annually. There are three service intervals: Annual, Intermediate and Major Servicing.

Annual Servicing (2,000 Hrs/1 Year)

- Suction inlet filter o-ring, Lubri-jet oil (for o-ring) and system check

Intermediate Servicing (4,000 Hrs/2 Years)

- Annual Service parts and Service Kit for Hydrocyclone Unit (Amalgam Separator)

Major Servicing (10,000 Hrs/5 Years)

- Annual Service parts and Service Kit for Suction Motor, Separator and Internal Hose Kit

ANNUAL SERVICING SCHEDULE

Micro SMART



Year	Minor Service Kits	Major Service Kits
1	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
2	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
3	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
4	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
5	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-MSMSK0 Micro SMART Major Service Kit Includes: (x1) C-025401-SK (Micro SMART Separator) (x1) C-014057-SK (Micro SMART Blower) (x1) C-MSIHSK (Internal Hose)
6	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
7	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
8	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
9	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
10	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	

ANNUAL SERVICING SCHEDULE

Micro SMART with Hydrocyclone



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	{x1} C-211956 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)		
2	{x1} C-211956 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)	{x1} C-MR60SK ISO 6 Hydrocyclone	
3	{x1} C-211956 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)		
4	{x1} C-211956 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)	{x1} C-MR60SK ISO 6 Hydrocyclone	
5	{x1} C-211956 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)		{x1} C-MSMSK0 Micro SMART Major Service Kit Includes: {x1} C-025401-SK (Micro SMART Separator) {x1} C-014057-SK (Micro SMART Blower) {x1} C-MSIHSK (Internal Hose)
6	{x1} C-211956 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)	{x1} C-MR60SK ISO 6 Hydrocyclone	
7	{x1} C-211956 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)		
8	{x1} C-211956 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)	{x1} C-MR60SK ISO 6 Hydrocyclone	
9	{x1} C-211956 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)		
10	{x1} C-211956 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)	{x1} C-MR60SK ISO 6 Hydrocyclone	

ANNUAL SERVICING SCHEDULE

Micro SMART Cube



Year	Minor Service Kits	Major Service Kits
1	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
2	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
3	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
4	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
5	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-MSCMSK Micro SMART CUBE Major Service Kit Includes: (x1) C-025401-SK (Micro SMART Separator) (x1) C-017500-SK (Micro SMART Cube Blower) (x1) C-MSCIH0-SK (Internal Hose)
6	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
7	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
8	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
9	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
10	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	

ANNUAL SERVICING SCHEDULE

Micro SMART Cube with Hydrocyclone



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
2	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-MR60SK ISO 6 Hydrocyclone	
3	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
4	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-MR60SK ISO 6 Hydrocyclone	
5	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		(x1) C-MSCMSK Micro SMART CUBE Major Service Kit Includes: (x1) C-025401-SK (Micro SMART Separator) (x1) C-017500-SK (Micro SMART Cube Blower) (x1) C-MSCIH0-SK (Internal Hose)
6	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-MR60SK ISO 6 Hydrocyclone	
7	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
8	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-MR60SK ISO 6 Hydrocyclone	
9	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
10	(x1) C-211956 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-MR60SK ISO 6 Hydrocyclone	

ANNUAL SERVICING SCHEDULE

Turbo SMART A



Year	Minor Service Kits	Major Service Kits
1	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
2	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
3	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
4	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
5	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-TSMSK0 Turbo SMART Major Service Kit Includes: (x1) C-0034MS-3SK (Turbo SMART Separator) (x1) C-1M0000-SK (Turbo SMART Blower) (x1) C-TSIHSK (Internal Hose)
6	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
7	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
8	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
9	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
10	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	

ANNUAL SERVICING SCHEDULE

Turbo SMART A with Hydrocyclone



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
2	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-1H00SK ISO18 Hydrocyclone	
3	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
4	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-1H00SK ISO18 Hydrocyclone	
5	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		(x1) C-TSMSK0 Turbo SMART Major Service Kit Includes: (x1) C-0034MS-3SK (Turbo SMART Separator) (x1) C-1M0000-SK (Turbo SMART Blower) (x1) C-TSIHSK (Internal Hose)
6	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-1H00SK ISO18 Hydrocyclone	
7	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
8	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-1H00SK ISO18 Hydrocyclone	
9	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
10	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(1x) C-1H00SK ISO18 Hydrocyclone	

ANNUAL SERVICING SCHEDULE

Turbo SMART B



Year	Minor Service Kits	Major Service Kits
1	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
2	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
3	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
4	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
5	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-TSMSK0 Turbo SMART Major Service Kit Includes: (x1) C-0034MS-3SK (Turbo SMART Separator) (x1) C-1M0000-SK (Turbo SMART Blower) (x1) C-TSIHSC (Internal Hose)
6	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
7	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
8	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
9	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
10	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	

ANNUAL SERVICING SCHEDULE


Turbo SMART B with Hydrocyclone



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	{x1} C-211974 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)		
2	{x1} C-211974 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)	{x1} C-1H00SK ISO18 Hydrocyclone	
3	{x1} C-211974 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)		
4	{x1} C-211974 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)	{x1} C-1H00SK ISO18 Hydrocyclone	
5	{x1} C-211974 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)		{x1} C-TSMK0 Turbo SMART Major Service Kit Includes: {x1} C-0034MS-3SK (Turbo SMART Separator) {x1} C-1M0000-SK (Turbo SMART Blower) {x1} C-TSIHSK (Internal Hose)
6	{x1} C-211974 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)	{x1} C-1H00SK ISO18 Hydrocyclone	
7	{x1} C-211974 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)		
8	{x1} C-211974 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)	{x1} C-1H00SK ISO18 Hydrocyclone	
9	{x1} C-211974 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)		
10	{x1} C-211974 Inlet filter o-ring {x1} M-040727 Lubrijet (single bottle)	{1x} C-1H00SK ISO18 Hydrocyclone	

ANNUAL SERVICING SCHEDULE

Turbo SMART A Cube

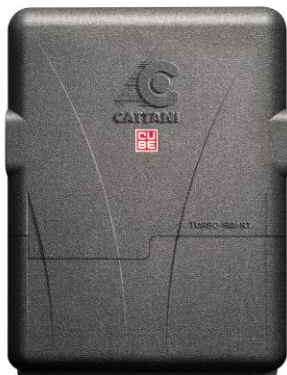

Year	Minor Service Kits	Major Service Kits
1	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
2	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
3	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
4	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
5	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-TSCMSK Turbo SMART Cube Major Service Kit Includes: (x1) C-0034MS-3SK (Turbo SMART Separator) (x1) C-014095-SK (Turbo SMART Cube Blower) (x1) C-TSCIH0-SK (Internal Hose)
6	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
7	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
8	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
9	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
10	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	

ANNUAL SERVICING SCHEDULE

Turbo SMART A Cube with Hydrocyclone



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
2	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-1H00SK ISO18 Hydrocyclone	
3	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
4	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-1H00SK ISO18 Hydrocyclone	
5	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		(x1) C-TSCMSK Turbo SMART Cube Major Service Kit Includes: (x1) C-0034MS-3SK (Turbo SMART Separator) (x1) C-014095-SK (Turbo SMART Cube Blower) (x1) C-TSCIH0-SK (Internal Hose)
6	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-1H00SK ISO18 Hydrocyclone	
7	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
8	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-1H00SK ISO18 Hydrocyclone	
9	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
10	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(1x) C-1H00SK ISO18 Hydrocyclone	

ANNUAL SERVICING SCHEDULE


Turbo SMART B Cube



Year	Minor Service Kits	Major Service Kits
1	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
2	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
3	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
4	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
5	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-TSCMSK Turbo SMART Cube Major Service Kit Includes: (x1) C-0034MS-3SK (Turbo SMART Separator) (x1) C-014095-SK (Turbo SMART Cube Blower) (x1) C-TSCIHO-SK (Internal Hose)
6	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
7	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
8	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
9	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	
10	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	

ANNUAL SERVICING SCHEDULE

Turbo SMART B Cube with Hydrocyclone

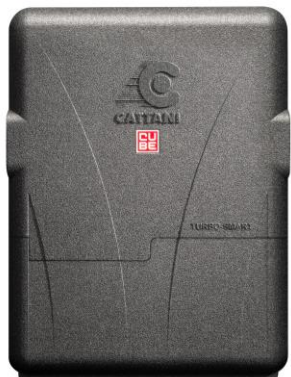

Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
2	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-1H00SK ISO18 Hydrocyclone	
3	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
4	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-1H00SK ISO18 Hydrocyclone	
5	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		(x1) C-TSCMSK Turbo SMART Cube Major Service Kit Includes: (x1) C-0034MS-3SK (Turbo SMART Separator) (x1) C-014095-SK (Turbo SMART Cube Blower) (x1) C-TSCIHO-SK (Internal Hose)
6	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-1H00SK ISO18 Hydrocyclone	
7	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
8	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-1H00SK ISO18 Hydrocyclone	
9	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
10	(x1) C-211974 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(1x) C-1H00SK ISO18 Hydrocyclone	

ANNUAL SERVICING SCHEDULE

Tandem Micro SMART



Year	Minor Service Kits	Major Service Kits
1	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
2	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
3	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
4	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
5	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-MSMSK0 Micro SMART Major Service Kit Includes: (x2) C-025401-SK (Micro SMART Separator) (x2) C-014057-SK (Micro SMART Blower) (x2) C-MSIHSK (Internal Hose)
6	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
7	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
8	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
9	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
10	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	

ANNUAL SERVICING SCHEDULE

Tandem Micro SMART with Hydrocyclone



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		
2	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-MR60SK- ISO 6 Hydrocyclone	
3	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		
4	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-MR60SK- ISO 6 Hydrocyclone	
5	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		(x2) C-MSMSK0 Micro SMART Major Service Kit Includes: (x2) C-025401-SK (Micro SMART Separator) (x2) C-014057-SK (Micro SMART Blower) (x2) C-MSIHSK (Internal Hose)
6	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-MR60SK- ISO 6 Hydrocyclone	
7	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		
8	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-MR60SK- ISO 6 Hydrocyclone	
9	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		
10	(x2) C-211956 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-MR60SK- ISO 6 Hydrocyclone	

ANNUAL SERVICING SCHEDULE

Tandem Turbo SMART A



Year	Minor Service Kits	Major Service Kits
1	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
2	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
3	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
4	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
5	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-TSMSK0 Turbo SMART Major Service Kit Includes: (x2) C-0034MS-3SK (Turbo SMART Separator) (x2) C-1M0000-SK (Turbo SMART Blower) (x2) C-TSIHSK (Internal Hose)
6	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
7	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
8	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
9	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
10	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	

ANNUAL SERVICING SCHEDULE

Tandem Turbo SMART A with Hydrocyclone



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		
2	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-1H00SK ISO18 Hydrocyclone	
3	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		
4	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-1H00SK ISO18 Hydrocyclone	
5	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		(x2) C-TSMSK0 Turbo SMART Major Service Kit Includes: (x2) C-0034MS-3SK (Turbo SMART Separator) (x2) C-1M0000-SK (Turbo SMART Blower) (x2) C-TSIHSK (Internal Hose)
6	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-1H00SK ISO18 Hydrocyclone	
7	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		
8	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-1H00SK ISO18 Hydrocyclone	
9	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		
10	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-1H00SK ISO18 Hydrocyclone	

ANNUAL SERVICING SCHEDULE

Tandem Turbo SMART B



Year	Minor Service Kits	Major Service Kits
1	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
2	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
3	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
4	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
5	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-TSMSK0 Turbo SMART Major Service Kit Includes: (x2) C-0034MS-3SK (Turbo SMART Separator) (x2) C-1M0000-SK (Turbo SMART Blower) (x2) C-TSIHSK (Internal Hose)
6	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
7	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
8	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
9	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	
10	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	

ANNUAL SERVICING SCHEDULE

Tandem Turbo SMART B with Hydrocyclone



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		
2	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-1H00SK ISO18 Hydrocyclone	
3	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		
4	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-1H00SK ISO18 Hydrocyclone	
5	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		(x2) C-TSMK0 Turbo SMART Major Service Kit Includes: (x2) C-0034MS-3SK (Turbo SMART Separator) (x2) C-1M0000-SK (Turbo SMART Blower) (x2) C-TSIHSK (Internal Hose)
6	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-1H00SK ISO18 Hydrocyclone	
7	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		
8	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-1H00SK ISO18 Hydrocyclone	
9	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)		
10	(x2) C-211974 Inlet filter o-ring (x2) M-040727 Lubrijet (single bottle)	(x2) C-1H00SK ISO18 Hydrocyclone	

ANNUAL SERVICING SCHEDULE

Multi SMART T3




Year	Minor Service Kits	Major Service Kits
1	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	
2	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	
3	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	
4	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	
5	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	(x3) C-TSMSK0 Turbo SMART Major Service Kit Includes: (x3) C-0034MS-3SK (Turbo SMART Separator) (x3) C-1M0000-SK (Turbo SMART Blower) (x3) C-TSIHSK (Internal Hose)
6	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	
7	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	
8	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	
9	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	
10	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	

* Replace (3x) C-160540 annually if HEPA Filters are fitted.

ANNUAL SERVICING SCHEDULE

Multi SMART T3 with Hydrocyclone

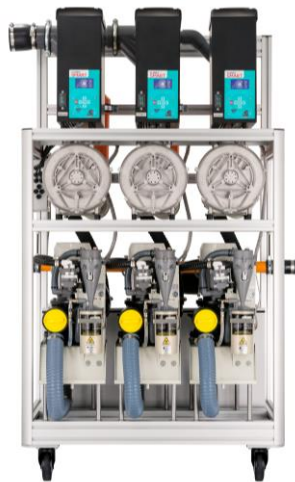

Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)		
2	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	(x3) C-1H00SK ISO18 Hydrocyclone	
3	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)		
4	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	(x3) C-1H00SK ISO18 Hydrocyclone	
5	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)		(x3) C-TSMSK0 Turbo SMART Major Service Kit Includes: (x3) C-0034MS-3SK (Turbo SMART Separator) (x3) C-1M0000-SK (Turbo SMART Blower) (x3) C-TSIHSK (Internal Hose)
6	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	(x3) C-1H00SK ISO18 Hydrocyclone	
7	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)		
8	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	(x3) C-1H00SK ISO18 Hydrocyclone	
9	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)		
10	(x3) C-211974 Inlet filter o-ring (x3) M-040727 Lubrijet (single bottle)	(x3) C-1H00SK ISO18 Hydrocyclone	

* Replace (3x) C-160540 annually if HEPA Filters are fitted.

ANNUAL SERVICING SCHEDULE

Multi SMART T4



Year	Minor Service Kits	Major Service Kits
1	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle)	
2	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle)	
3	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle)	
4	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle)	
5	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle)	(x4) C-TSMK0 Turbo SMART Major Service Kit Includes: (x4) C-0034MS-3SK (Turbo SMART Separator) (x4) C-1M0000-SK (Turbo SMART Blower) (x4) C-TSIHSK (Internal Hose)
6	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle)	
7	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle)	
8	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle)	
9	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle)	
10	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle)	

* Replace (4x) C-160540 annually if HEPA Filters are fitted.

ANNUAL SERVICING SCHEDULE

Multi SMART T4 with Hydrocyclone



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	{x4} C-211974 Inlet filter o-ring {x4} M-040727 Lubrijet (single bottle)		
2	{x4} C-211974 Inlet filter o-ring {x4} M-040727 Lubrijet (single bottle)	{x4} C-1H00SK ISO18 Hydrocyclone	
3	{x4} C-211974 Inlet filter o-ring {x4} M-040727 Lubrijet (single bottle)		
4	{x4} C-211974 Inlet filter o-ring {x4} M-040727 Lubrijet (single bottle)	{x4} C-1H00SK ISO18 Hydrocyclone	
5	{x4} C-211974 Inlet filter o-ring {x4} M-040727 Lubrijet (single bottle)		{x4} C-TSMSK0 Turbo SMART Major Service Kit Includes: {x4} C-0034MS-3SK (Turbo SMART Separator) {x4} C-1M0000-SK (Turbo SMART Blower) {x4} C-TSIHSK (Internal Hose)
6	{x4} C-211974 Inlet filter o-ring {x4} M-040727 Lubrijet (single bottle)	{x4} C-1H00SK ISO18 Hydrocyclone	
7	{x4} C-211974 Inlet filter o-ring {x4} M-040727 Lubrijet (single bottle)		
8	{x4} C-211974 Inlet filter o-ring {x4} M-040727 Lubrijet (single bottle)	{x4} C-1H00SK ISO18 Hydrocyclone	
9	{x4} C-211974 Inlet filter o-ring {x4} M-040727 Lubrijet (single bottle)		
10	{x4} C-211974 Inlet filter o-ring {x4} M-040727 Lubrijet (single bottle)	{x4} C-1H00SK ISO18 Hydrocyclone	

* Replace (4x) C-160540 annually if HEPA Filters are fitted.

ANNUAL SERVICING SCHEDULE

Multi SMART T5



Year	Minor Service Kits	Major Service Kits
1	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	
2	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	
3	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	
4	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	
5	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	(x5) C-TSMSK0 Turbo SMART Major Service Kit Includes: (x5) C-0034MS-3SK (Turbo SMART Separator) (x5) C-1M0000-SK (Turbo SMART Blower) (x5) C-TSIHSK (Internal Hose)
6	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	
7	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	
8	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	
9	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	
10	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	

* Replace (5x) C-160540 annually if HEPA Filters are fitted.

ANNUAL SERVICING SCHEDULE

Multi SMART T5 with Hydrocyclone



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)		
2	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	(x5) C-1H00SK ISO18 Hydrocyclone	
3	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)		
4	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	(x5) C-1H00SK ISO18 Hydrocyclone	
5	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)		(x5) C-TSMSK0 Turbo SMART Major Service Kit Includes: (x5) C-0034MS-3SK (Turbo SMART Separator) (x5) C-1M0000-SK (Turbo SMART Blower) (x5) C-TSIHSK (Internal Hose)
6	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	(x5) C-1H00SK ISO18 Hydrocyclone	
7	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)		
8	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	(x5) C-1H00SK ISO18 Hydrocyclone	
9	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)		
10	(x5) C-211974 Inlet filter o-ring (x5) M-040727 Lubrijet (single bottle)	(x5) C-1H00SK ISO18 Hydrocyclone	

* Replace (5x) C-160540 annually if HEPA Filters are fitted.

ANNUAL SERVICING SCHEDULE

**Maxi SMART
with Hydrocyclone
and Acoustic Cover**



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
2	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-043390-SK (ISO60 Hydrocyclone)	
3	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
4	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-043390-SK (ISO60 Hydrocyclone)	
5	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		(x1) C-MAXIMSK Maxi SMART Major Service Kit (x1) C-025408-SK (Maxi SMART Separator) (x1) C-091655-SK (Maxi SMART Blower) (x1) C-MXIHSK (Internal Hose)
6	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-043390-SK (ISO60 Hydrocyclone)	
7	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
8	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-043390-SK (ISO60 Hydrocyclone)	
9	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
10	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-043390-SK (ISO60 Hydrocyclone)	

ANNUAL SERVICING SCHEDULE

Maxi SMART with
hydrocyclone
without Acoustic
Cover



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
1	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
2	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-043390-SK (ISO60 Hydrocyclone)	
3	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
4	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-043390-SK (ISO60 Hydrocyclone)	
5	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		(x1) C-MAXIMSK Maxi SMART Major Service Kit Includes: (x1) C-025408-SK (Maxi SMART Separator) (x1) C-091655-SK (Maxi SMART Blower) (x1) C-MXIHSK (Internal Hose)
6	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-043390-SK (ISO60 Hydrocyclone)	
7	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
8	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-043390-SK (ISO60 Hydrocyclone)	
9	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)		
10	(x1) C-215986 Inlet filter o-ring (x1) M-040727 Lubrijet (single bottle)	(x1) C-043390-SK (ISO60 Hydrocyclone)	

ANNUAL SERVICING SCHEDULE

Blok Jet Large Central Vacuum System (200 Chairs)



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
0.5	(x2) C-AS1000 - Antifoam Liquid		
1	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle) (x2) C-AS1000 - Antifoam Liquid (x3) C-042094 -HEPA H14 Filter Cartridge		
1.5	(x2) C-AS1000 -Antifoam Liquid		
2	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle) (x2) C-AS1000 - Antifoam Liquid (x3) C-042094 -HEPA H14 Filter Cartridge	(x4) C-1H00SK ISO18 Hydrocyclone	
2.5	(x2) C-AS1000 - Antifoam Liquid		
3	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle) (x2) C-AS1000 - Antifoam Liquid (x3) C-042094 -HEPA H14 Filter Cartridge		
3.5	(x2) C-AS1000 - Antifoam Liquid		
4	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle) (x2) C-AS1000 - Antifoam Liquid (x3) C-042094 -HEPA H14 Filter Cartridge	(x4) C-1H00SK ISO18 Hydrocyclone	
4.5	(x2) C-AS1000 - Antifoam Liquid		
5	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle) (x2) C-AS1000 - Antifoam Liquid (x3) C-042094 -HEPA H14 Filter Cartridge		(x4) C-083067-SK (Suction Pump – Blokjet) (x4) C-TSMSK0 Turbo SMART Major Service Kit Includes: (x4) C-0034MS-3SK (Turbo SMART Separator) (x4) C-1M0000-SK (Turbo SMART Blower) (x4) C-TSIHSK (Internal Hose)

ANNUAL SERVICING SCHEDULE

Blok Jet Large Central Vacuum System (200 Chairs)



Year	Minor Service Kits	Intermediate Service Kits	Major Service Kits
5.5	(x2) C-AS1000 - Antifoam Liquid		
6	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle) (x2) C-AS1000 - Antifoam Liquid (x3) C-042094 -HEPA H14 Filter Cartridge		
6.5	(x2) C-AS1000 -Antifoam Liquid		
7	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle) (x2) C-AS1000 - Antifoam Liquid (x3) C-042094 -HEPA H14 Filter Cartridge	(x4) C-1H00SK ISO18 Hydrocyclone	
7.5	(x2) C-AS1000 - Antifoam Liquid		
8	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle) (x2) C-AS1000 - Antifoam Liquid (x3) C-042094 -HEPA H14 Filter Cartridge		
8.5	(x2) C-AS1000 - Antifoam Liquid		
9	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle) (x2) C-AS1000 - Antifoam Liquid (x3) C-042094 -HEPA H14 Filter Cartridge	(x4) C-1H00SK ISO18 Hydrocyclone	
9.5	(x2) C-AS1000 - Antifoam Liquid		
10	(x4) C-211974 Inlet filter o-ring (x4) M-040727 Lubrijet (single bottle) (x2) C-AS1000 - Antifoam Liquid (x3) C-042094 -HEPA H14 Filter Cartridge		(x4) C-083067-SK (Suction Pump – Blokjet) *Assess Turbo SMART for life expectancy

ANNUAL SERVICING SCHEDULE

Aspi Jet 6



Year	Minor Service Kits	Major Service Kits
1	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
2	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
3	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
4	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
5	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	(x1) C-002300-SK (Blower Uni Jet 75 ~1)
6	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
7	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
8	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
9	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
10	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	

ANNUAL SERVICING SCHEDULE

Aspi Jet 6

Fitted with Bacterial
Filter



Year	Minor Service Kits	Major Service Kits
1	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
2	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
3	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
4	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
5	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	(x1) C-002300-SK (Blower Uni Jet 75 ~1)
6	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
7	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
8	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
9	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
10	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	

ANNUAL SERVICING SCHEDULE

Aspi Jet 7



Year	Minor Service Kits	Major Service Kits
1	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
2	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
3	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
4	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
5	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	(x1) C-002300-SK (Blower Uni Jet 75 ~1)
6	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
7	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
8	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
9	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	
10	(x1) C-000600-SK (x1) M-040727 Lubrijet (single bottle)	

ANNUAL SERVICING SCHEDULE

Aspi Jet 7
Fitted with
Bacterial Filter



Year	Minor Service Kits	Major Service Kits
1	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
2	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
3	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
4	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
5	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	(x1) C-002300-SK (Blower Uni Jet 75 ~1)
6	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
7	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
8	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
9	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	
10	(x1) C-000600-SKB (x1) M-040727 Lubrijet (single bottle)	

ANNUAL SERVICING SCHEDULE

Turbo Jet 1 Modular



Year	Major Service Kits
1	
2	
3	
4	
5	{x1} C-002300-SK (Blower Uni Jet 75 ~1) {x1} C-0033M0-SK (T1/T2 Separator ~1)
6	
7	
8	
9	
10	

ANNUAL SERVICING SCHEDULE

Turbo Jet 2 Modular



Year	Major Service Kits
1	
2	
3	
4	
5	(x1) C-0023V0-SK (Blower Uni Jet 75 2V ~1) (x1) C-0033M0-SK (T1/T2 Separator ~1)
6	
7	
8	
9	
10	

ANNUAL SERVICING SCHEDULE

Uni Jet 75 - 0.4 kW

Single-Phase



Year	Major Service Kits
1	
2	
3	
4	
5	(x1) C-002300-SK (Blower Uni Jet 75 ~1)
6	
7	
8	
9	
10	

ANNUAL SERVICING SCHEDULE

Uni Jet 75 2V - 0.7 kW

Single-Phase



Year	Major Service Kits
1	
2	
3	
4	
5	(x1) C-0023V0-SK (Blower Uni Jet 75 2V ~1)
6	
7	
8	
9	
10	

ANNUAL SERVICING SCHEDULE

Tecno Jet - 0.75 kW
Single-Phase



Year	Major Service Kits
1	
2	
3	
4	
5	(x1) C-006900-SK (Blower Tecno Jet ~1)
6	
7	
8	
9	
10	

ANNUAL SERVICING SCHEDULE

Flux Jet - 1.1 kW
Single-Phase



Year	Major Service Kits
1	
2	
3	
4	
5	(x1) C-014000-SK (Blower Flux Jet ~1)
6	
7	
8	
9	
10	

ANNUAL SERVICING SCHEDULE

Medio Jet - 3 kW

Three-Phase



Year	Major Service Kits
1	
2	
3	
4	
5	(x1) C-A00150-SK (Blower Medio Jet ~3)
6	
7	
8	
9	
10	

**DISCLAIMER**

Failure to observe a danger sign or warning may harm operators, patients or technical staff. Warning signs cannot always fully express all hazards and dangers associated with equipment repairs.

Prior to commencing work, ensure all appropriate Personal Protective Equipment (PPE) such as: gloves, mask, full sleeve gown and eye protection are worn. Refer to the product "User Manual" for manufacturer's recommendations.

Cattani Technical Support:

+61 3 9484 1120

cattani@cattani.com.au



SALES AND TECHNICAL SUPPORT:

+61 3 9484 1120
CATTANI@CATTANI.COM.AU

COMPANY ADDRESS:

280 DUNDAS STREET
THORBURY VIC 3071

WWW.CATTANI.COM.AU