


RAISING THE BAR

Personalised Premium Dental Suction

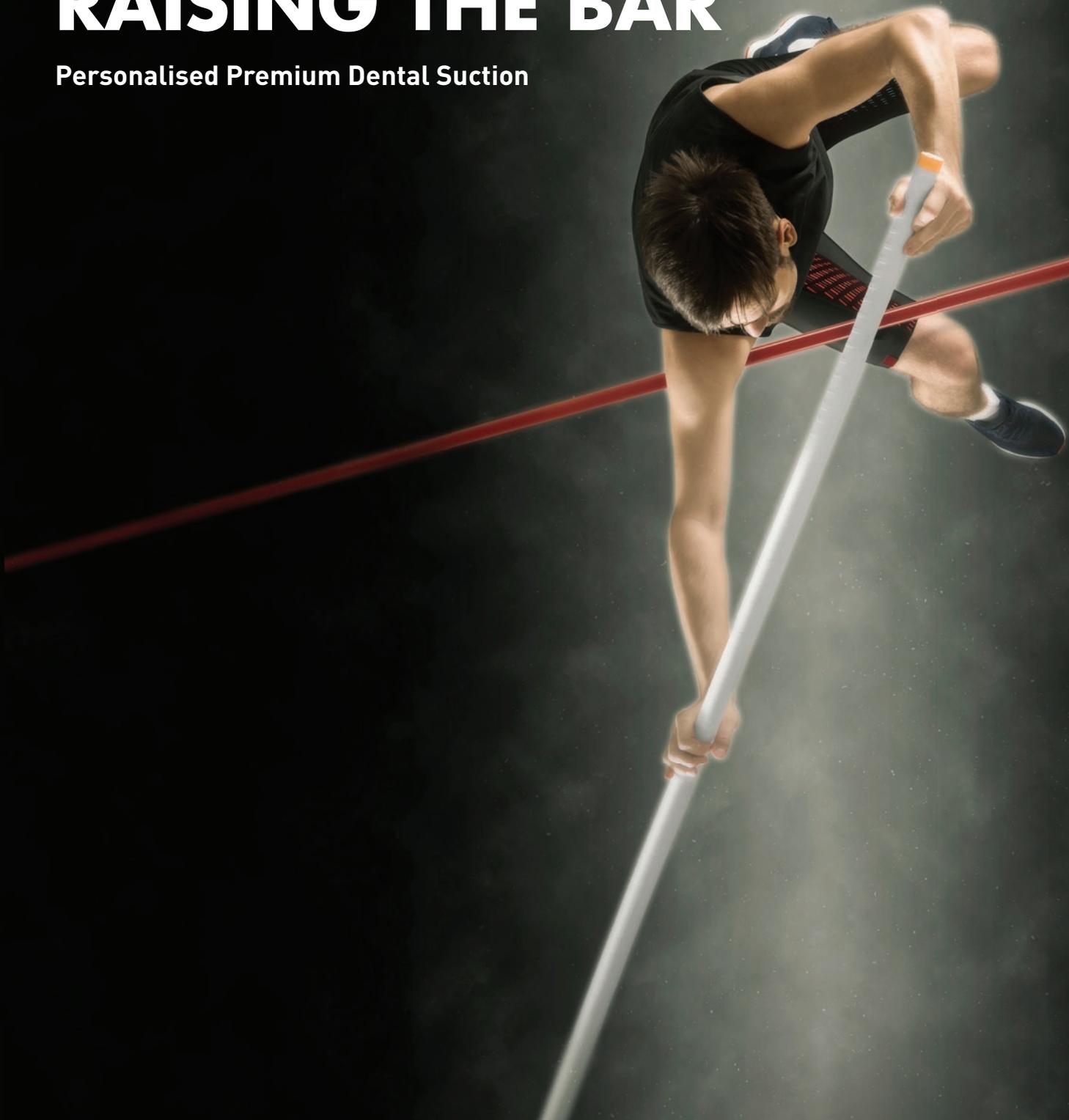

NEW TURBO SMART HP

High performance. Highly personalised.

From dental hygiene to oral surgery, the all-new TURBO SMART HP is the next generation of first-class dental suction to deliver constant peak performance, fully customisable to suit your dental clinic and individual requirements.



NEW TURBO SMART HP

High-Performance | 338 mbar Vacuum Pressure | 3 Chairs*



Turbo SMART HP with ISO 18 Hydrocyclone

Item Code 035145

SPECIFICATIONS

Power 2.7kW 12A
1PH 230V

Sound Output dB(A) 70

Max Air Flow at l/min 1916

Max Vacuum for Continuous Operation mBar 338

DIMENSIONS

Height mm 765

Width mm 460

Depth mm 350

Weight kg 55

Details That Matter

- High performance 338 mBar vacuum pressure for highly personalised dental suction performance
- Capability to customise suction power to serve 4 surgeries*
- Purpose-built for intraoral dental aerosol control
- 70dB(A) – improved noise reduction via ventilation technology
- Improved dual-contour blade design for more suction power while using less energy

Compare SMART Series



	Micro SMART	Turbo SMART A	Turbo SMART B	Turbo SMART HP
No. of Chairs	2 Chairs	3 Chairs	4 Chairs	3 Chairs*
Maximum Flow	900 l/min	1400 l/min	1750 l/min	1916 l/min
Maximum Head for Continuous Duty	210 mBar	210 mBar	210 mBar	338 mBar
Sound Pressure Without Box	64 dB(A)	69 dB(A)	73.7 dB(A)	70dB(A)
Amalgam Separator	Optional	Optional	Optional	Standard
Suction Motor	Single (120Hz)	Single (85Hz)	Single (110Hz)	Twin in Series (90Hz)
Dual Inverter	Yes	Yes	Yes	Yes
Suction Impeller	Original SMART	Original SMART	Original SMART	Patented SMART dual-contour blade
Quick-Connect Cable Connection	No	No	No	Standard
Ventilation Cooling	3 Fans	3 Fans	3 Fans	4 Fans

* 3 chairs as standard, 4 chairs dependent on maximum head setting. For more information, please contact us on +61 3 9484 1120.

Sales and Technical Support



Australia – Phone +61 3 9484 1120 Email cattani@cattani.com.au Web www.cattani.com.au

New Zealand – Phone 0800 68 22 88 Email cattani@cattani.co.nz Web www.cattani.co.nz

Company Address: 280 Dundas Street, Thornbury VIC 3071